


clinitherm[®] pasteur

Pasteurization of breast milk
– safe and reliable

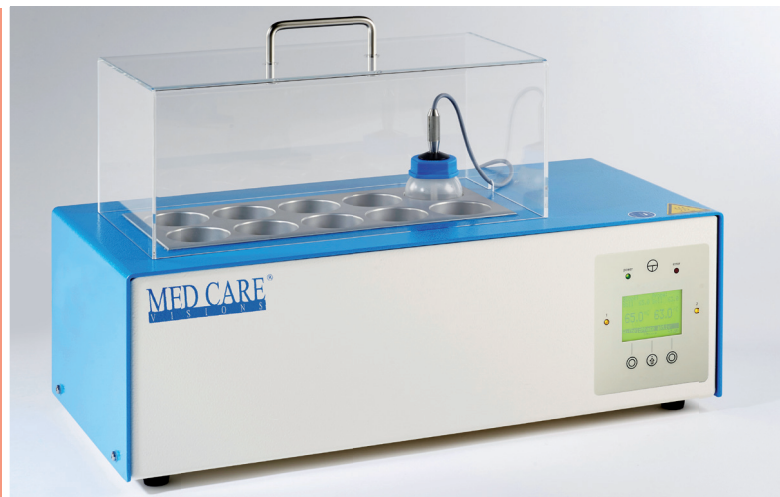

The clinitherm[®] pasteur is a dry warming device used for pasteurizing breast milk.

Pasteurization

The pasteurization parameters in the clinitherm[®] pasteur are set at the manufacturer. Breast milk is pasteurized at a temperature of +63°C for 30 minutes. A reference milk bottle measures the temperature.

Safety

- Dry warming
- Visual and acoustic signal after the pasteurization process
- Clear and unambiguous visualization due to clear insulation cover
- Documentation software (optional)

Benefits:

- Germ free pasteurization due to dry warming
- Threefold safety due to software-based, electronic and mechanical overtemperature protection
- Permanent temperature control
- Menu-controlled user guidance
- Full automatic processing
- No use of disposables – saves costs
- Can also be used as a baby bottle warmer
- Available for 6 and 10 bottles



Reference bottle



Sensor

Ordering Information

Product Description	Ord.-No.
Clinitherm® pasteur, 6 bottles, 80mm depth, 54mm diameter	BK21410084
clinitherm® pasteur, 6 bottles, 80mm depth, 56mm diameter	BK21410084-56
Clinitherm® pasteur, 6 bottles, 120mm depth, 54mm diameter	BK21410092
Clinitherm® pasteur, 6 bottles, 120mm depth, 56mm diameter	BK21410092-56
Clinitherm® pasteur, 10 bottles, 80mm depth, 54mm diameter	BK21410086
Clinitherm® pasteur, 10 bottles, 80mm depth, 56mm diameter	BK21410086-56
Clinitherm® pasteur, 10 bottles, 120mm depth, 54mm diameter	BK21410094
Clinitherm® pasteur, 10 bottles, 120mm depth, 56mm diameter	BK21410094-56
Passive cooling block / Equipment	
Passive cooling block, for 6 bottles, Ø ca. 54mm depth 80mm	BK21410138
Passive cooling block for 6 bottles, Ø ca. 54mm depth 120mm	BK21410139
Passive cooling block for 6 bottles, Ø ca. 56mm depth 120mm	BK21410140
Passive cooling block for 10 bottles, Ø ca. 54mm depth 80mm	BK21410141
Passive cooling block for 10 bottles, Ø ca. 54mm depth 120mm	BK21410142
Passive cooling block for 10 bottles, Ø ca. 56mm depth 120mm	BK21410143
Thermostat brush for cleaning	BK21510040
Software clinitherm® pasteur	BK21410075
Reference probe for clinitherm®	BK21410125
Clinitherm® Baby/Pasteur accessory set (1x thermostat brush, 1x PU Dry-Pads Clinitherm®, 1x Pads gripper)	BK21510042